Biatain® Ibu

Indications:

Biatain® Ibu foam dressings are indicated for a wide range of exudating wounds, such as leg ulcers, pressure ulcers, diabetic foot ulcers, second degree burns, donor sites, postoperative wounds, and skin abrasions.

Biatain® Ibu foam dressings provide moist wound healing and may reduce wound pain caused by tissue damage.

Biatain® Ibu foam dressings are suitable for use in combination with compression therapy.

Directions for Use:

1. Rinse the wound according to clinical practice. To optimize the use of the dressing, select a dressing that overlaps the wound edge by a minimum of two centimeters.
2. **Biatain® Ibu** must be applied with the plain (non-printed) side to the wound.
3. **Biatain® Ibu** should be changed when clinically indicated or when visible signs of exudate approach the edge of the dressing. Depending upon the condition of the wound the dressing may be left in place for up to seven days.
4. The semi-permeable top film is shower proof and provides a barrier function to bacteria.
5. If the wound appears to show signs of deterioration please consult your Wound Care Specialist.

Biatain® Ibu


Non-Adhesive Foam Dressings with Ibuprofen

	size (cm)	item no.	available
	10x10	4110	☐
	10x20	4112	☐
	15x15	4115	☐
	20x20	4120	☐

Soft-Hold Foam Dressings with Ibuprofen

	size (cm)	item no.	available
	10x10	4140	☐
	10x20	4142	☐

Biatain® Ibu foam dressings combine superior exudate management and continuous release of ibuprofen^(1,2) locally with no observed systemic effect⁽²⁾. One dressing (10 x 10 cm) contains 50 mg of ibuprofen which is continuously released for up to 7 days in proportion to the level of exudate.

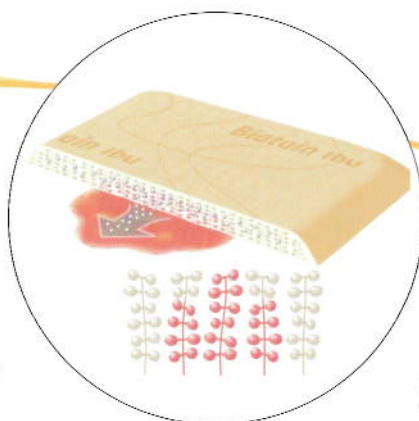Biatain

Coloplast Canada Corporation
3300 Ridgeway Dr. Unit 12 Mississauga, ON L5L 5Z9
www.coloplast.ca 1 888-880-8605

Coloplast and Biatain are registered trademarks of Coloplast A/S and Coloplast Canada Corporation
© 2009-06. All rights reserved Coloplast A/S, 3650 Hømløbsvej, Denmark

References

1. Jorgensen, B.; Friis, G. J.; Gottrup, F. Pain and quality of life for patients with venous leg ulcers: Proof of concept of the efficacy of Biatain Ibu, a new pain reducing wound dressing. Wound repair and regeneration, 2006, 14 [3], 233-239.
2. Steffensen, B.; Herping, S. P. K. Novel wound models for characterizing the effects of exudates levels on the controlled release of ibuprofen from foam dressings. European Wound Management Association, Prague, Czech Republic, 2006.

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Références

1. Jorgensen, B.; Friis, G. J.; Gottrup, F. Pain and quality of life for patients with venous leg ulcers: Proof of concept of the efficacy of Biatain Ibu, a new pain reducing wound dressing. Wound repair and regeneration, 2006, 14(3), 233-239.
2. Steffensen, B.; Herping, S. P. K. Novel wound models for characterizing the effects of exudates levels on the controlled release of ibuprofen from foam dressings, European Wound Management Association, Prague, Czech Republic, 2006.